


Virtual slide tray panel

The Tray panel is a virtual store of slides, user can send slides to the Tray and keep them in Tray for the current session. User can open select slides from tray and open them in viewer tabs or in a grid or remove them from the Tray.

Tray section is part of the “Home” tab and consists of 4 buttons.



Send slide to tray

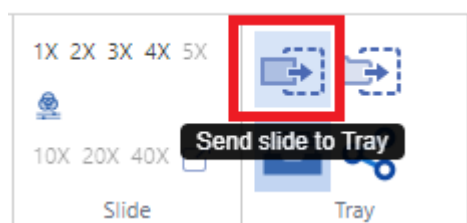


This button send opens slide to the tray panel. If a slide is already part of the tray it shows you a message that selected slide is already part of the tray.


Send directory to tray



This slide sends the whole selected directory of slides to tray. When user clicks the button a dialog box appears that lists out all the slides under that directory. It provides you detail about the selected server, directory itself, number of slides in it and current status of the tray items. All the slides appear in a checklist format with all slides selected by default, user can select/unselect individual slides as per need. The dialog box consists of couple buttons that includes button for appending selected slides to the tray panel. Other button “clear tray and send” clears the existing tray and adds selected slides to it.

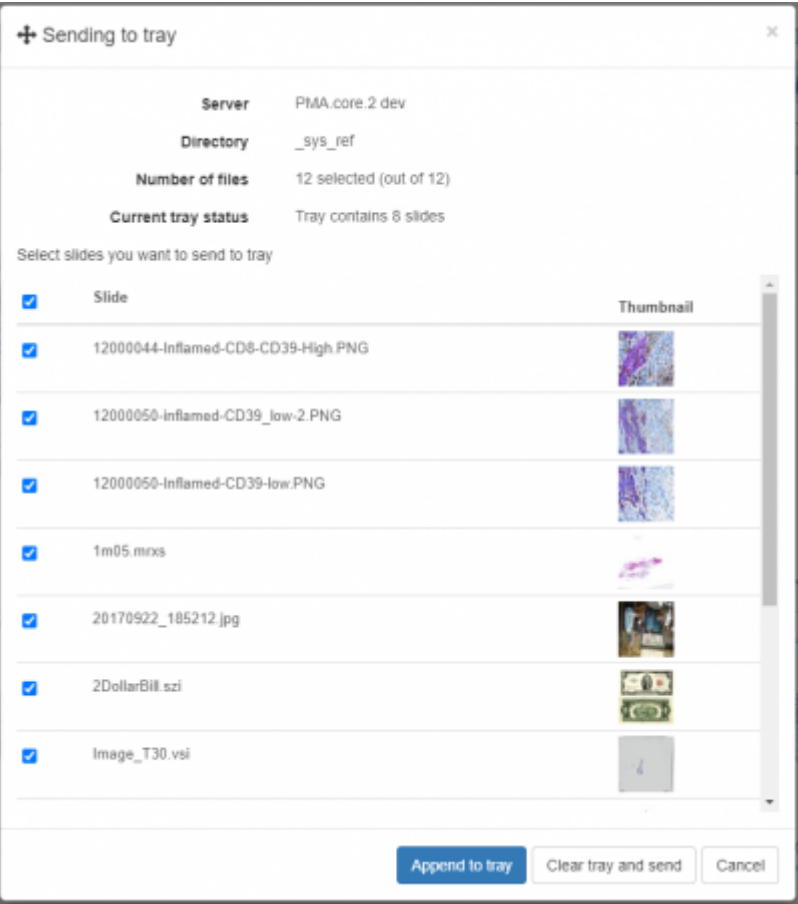

Here is how the slides appear in the tray panel:

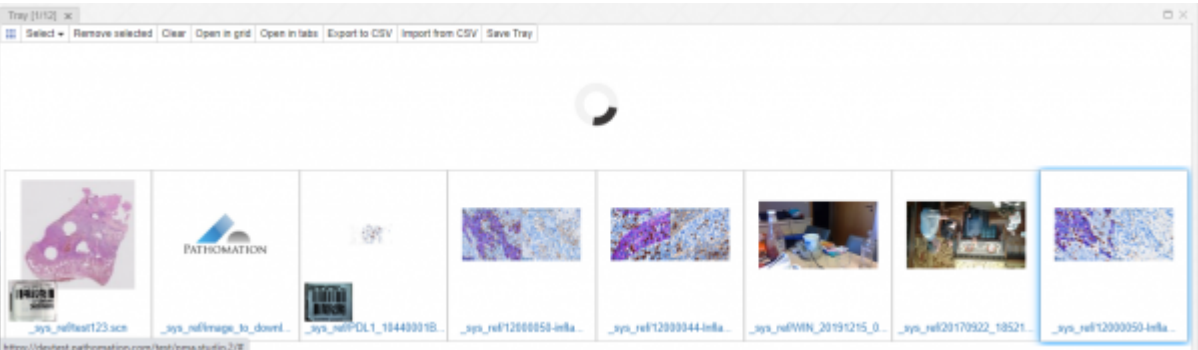

Tray panel layout

Tray panel can present slides in 2 ways - Meta data view and Slides view.


Metadata view

Metadata view displays slide information and shows all filled forms for the slide with it's values.


Slide info													
	Server	File name	Width	Height	Max Magnification	μ Per Pixel X	μ Per Pixel Y	Physical Size	Tile Size	Barcode	Overview	Layers	Channels
☐	PMA core.2 dev	_sys_refTest123.scn	36832	38432	20X	0.5	0.5	278.3 MB	512			1	1
☐	PMA core.2 dev	_sys_refImage_to_download.jpg	698	267	0	0	0	0.8 MB	512	N/A		1	1
☐	PMA core.2 dev	_sys_refPDL1_18440981860925.mrxs	88576	188624	48X	0.242797397769517	0.242797397769517	2491.6 MB	512			1	1
☐	PMA core.2 dev	_sys_ref12800058-Inflamed-CD38_low-2.PNG	1878	832	0	0	0	4.7 MB	512	N/A		1	1
☐	PMA core.2 dev	_sys_ref12000044-Inflamed-CD6-CD38-High.PNG	1880	833	0	0	0	4.7 MB	512	N/A		1	1
☐	PMA core.2 dev	_sys_refWIN_20191215_08_46_53_Pro.jpg	3264	1836	0	0	0	2.2 MB	512	N/A		1	1
☐	PMA core.2 dev	_sys_ref20170922_185212.jpg	2560	1440	0	0	0	1.3 MB	512	N/A		1	1

Slides view

Slides view shows slide contents in Gallery.

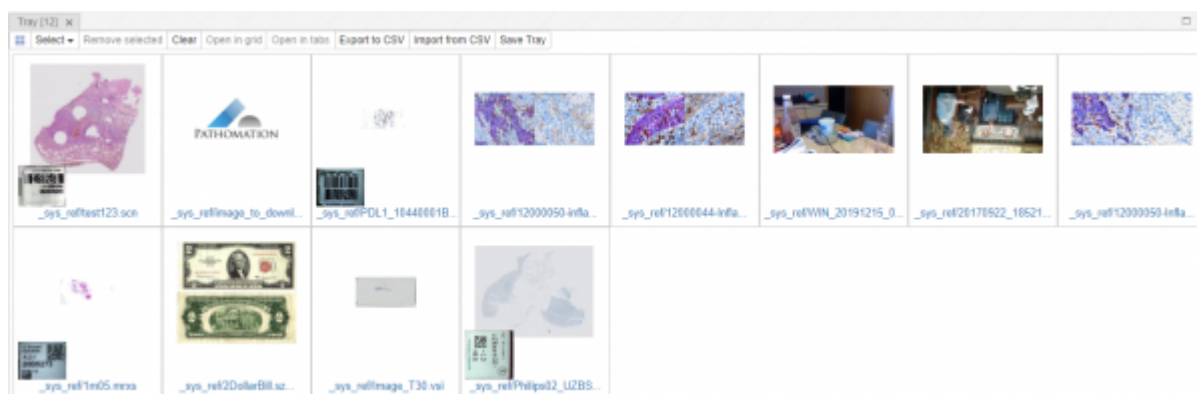

To make Tray persistent it's possible to export and import slides to/from tray in CSV format. Button toolbar is mostly the same for slides and metadata view except of 'Export to CSV' button which is available only for Metadata view. In both modes it's available to select slides by clicking on them, selected slides are being highlighted with blue color, in addition slide that is currently active in viewer will be additionally highlighted in Metadata view with blue dashed border.

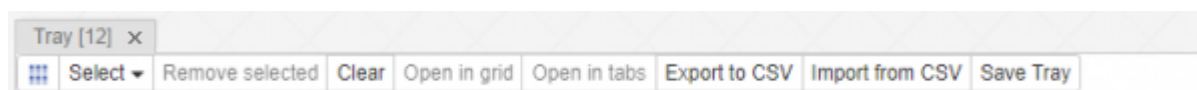

The Tray toolbar consist of several buttons, from the left to the right:

1. **Switch mode** - switches between Metadata and Slides mode
2. **Select dropdown** -
 1. All - select all slides in Tray
 2. None - deselect all slides in Tray - selection will be empty
 3. Invert - invert selection, deselect selected slides and vice versa
3. **Remove selected** - removes selected slides from the Tray
4. **Clear** - removes all slides from the Tray
5. **Open in grid** - open selected slides in one viewer panel in Grid, it's possible to perform this action only when selected slides count is less or equal 36 (6*6 grid)
6. **Open in tabs** - open each selected slide in separate viewer panel
7. **Export to CSV** - exports slides metadata information - slide info and filled forms to the csv file, download should start immediately

8. **Import from CSV** - opens an import wizard where you can drag an import file on it and it will load the valid slides to tray panel
9. **Save Tray** - saves all the changes made to tray and shows you a message if changes saved successfully. This allows you to save and then load the tray when you login again to the system

From:

<https://docs.pathomation.com/pma.studio/2.1.0/> - **PMA.studio 2.1.0**

Permanent link:

https://docs.pathomation.com/pma.studio/2.1.0/doku.php?id=panel_tray

Last update: **2023/04/06 21:52**

