



How to share a case?

1. Go to your cases via “My cases” on the sidebar.

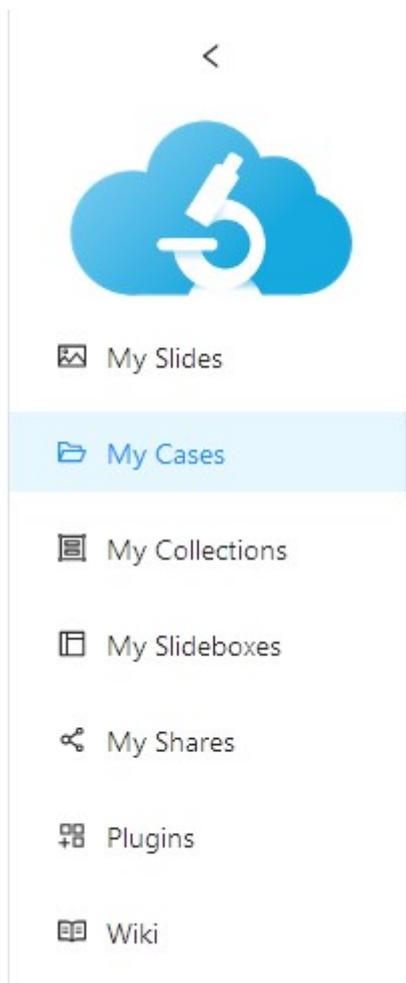

2. All your cases are shown in the gallery. On every case, you can see the share button. Click share button and a new window pops up.

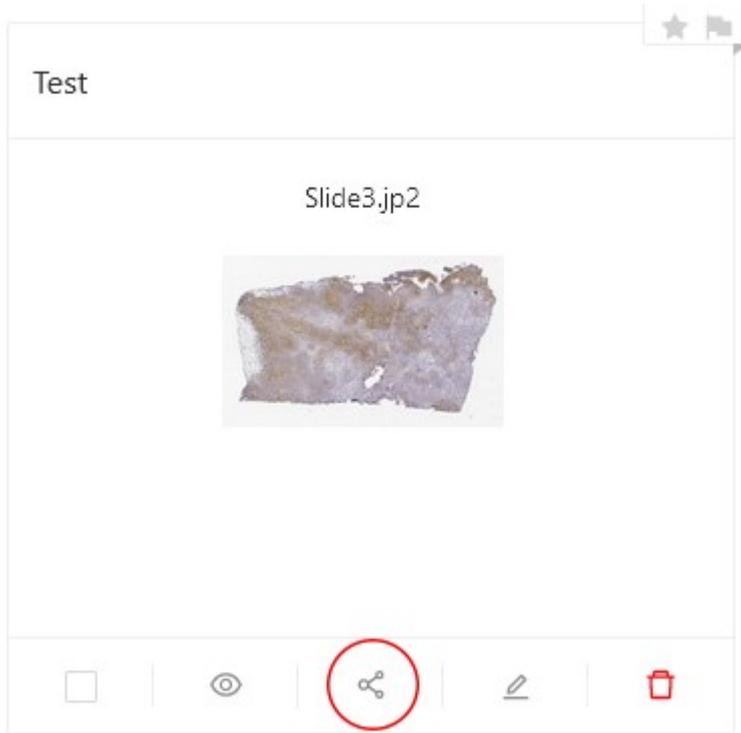

3. There are different ways to share your case.

1. You can copy the link and share it with others
2. You can embed your slide in your own web content
3. You can generate a QR-code for your slide
4. You can email the link of shared slide
5. You can directly send the link of shared slide via various social network service such as *Facebook, Twitter, WhatsApp, Messenger or LinkedIn.*

Also you can change the visibility option of your share. **Unlisted shares** can only be accessed by users with a valid link. **Public shares** can also be accessed by link, but are accessible from public listings too and anyone can view them.

Choose the options that suits your need best.

Share Test





Test
Created on Feb 10, 2022, 16:33
Contains 1 slide(s)

Share method

Visibility ⓘ: Public Unlisted

Link Embed QRcode Email Facebook Twitter WhatsApp Messenger

[ADVANCED \(more options and properties\)](#)

[Close](#)

4. By clicking on “ADVANCED (more options and properties)” you can finetune your share. You can password-protect your share, set an expiration date for your share or modify viewport to your needs.

Advanced sharing



Test
Created on Feb 10, 2022, 16:33
Contains 1 slide(s)

Share method

Visibility ⓘ: Public Unlisted

Link Embed QRcode Email Facebook Twitter WhatsApp Messenger LinkedIn

Share options


Password-protect the share:

Put an expiration on the share:

Viewport options



- Overview ⓘ: Visible Collapsed
- Dimensions control ⓘ: Visible Collapsed
- Barcode ⓘ: Visible Collapsed
- Rotation control ⓘ: Visible Collapsed
- Filename ⓘ: Visible Collapsed
- Scale line ⓘ: Visible Collapsed
- Snapshot button ⓘ: Visible
- Color adjustments ⓘ: Visible

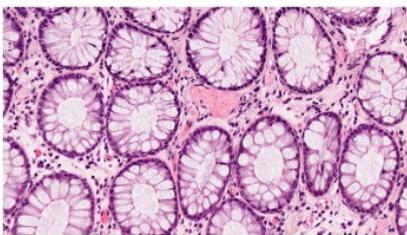

Channels
Default
Default
Time frame: 1 / 2
Z-stack: 2 / 2

Last update: 2022/02/14
15:22

how_to_share_a_case https://docs.pathomation.com/my/doku.php?id=how_to_share_a_case&rev=1644841353

From:
<https://docs.pathomation.com/my/> - **My Pathomation**

Permanent link:
https://docs.pathomation.com/my/doku.php?id=how_to_share_a_case&rev=1644841353

Last update: **2022/02/14 15:22**

