


# The viewer

PMA.control offers an efficient digital slide viewer that maximizes performance and ease of use. Whether you are using a desktop workstation with a computer mouse or a touch enabled device, you will find out that it is straight forward and smooth to use.

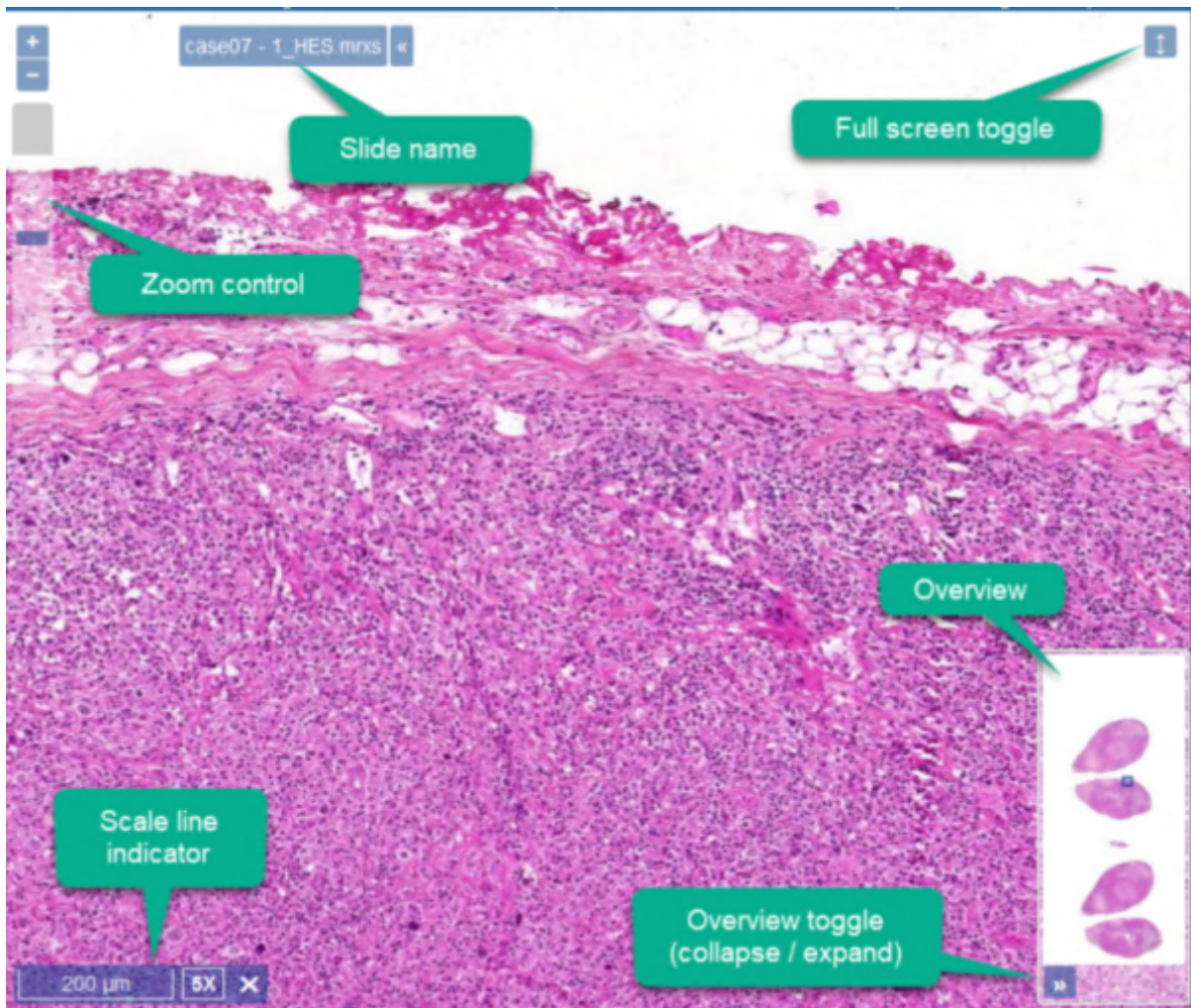

## Navigation

### MOUSE NAVIGATION

Pma.control works with all types of mice, however it is suggested that you use a mouse with a scroll wheel for the optimum performance. A mouse wheel with click-to-click scrolling is preferred for best user experience. The following diagram identifies how the mouse buttons are used.



## KEYBOARD + MOUSE

\* ALT + Mouse: Rotate \* SHIFT + Mouse: drag zoom

## KEYBOARD

- Press the LEFT ARROW, RIGHT ARROW, UP ARROW and DOWN ARROW keys.

## TOUCH ENABLED DEVICE

- zoom in and out by pinching the screen using two fingers

\* *Panning: To move the image around, drag it with the mouse or your finger if you are using a touch enabled device. You can also navigate by clicking within the overview image as well as by dragging the rectangle inside it. Horizontal and vertical panning Syncing viewers As already mentioned, it is possible to open multiple slides simultaneously. On top of this, you can “lock” the slides to each other, so that when you pan or zoom in one, the rest follow. This can be enabled by ticking the “Sync” check box, which can be found above the list of slides on the right slide of the page.*

From: <https://docs.pathomation.com/pathotrainer/> - Pathotrainer

Permanent link: [https://docs.pathomation.com/pathotrainer/doku.php?id=the\\_viewer&rev=1596811865](https://docs.pathomation.com/pathotrainer/doku.php?id=the_viewer&rev=1596811865)

Last update: 2020/08/07 17:51

