


Server settings

In this section, various aspects of the server behavior can be configured and controlled.

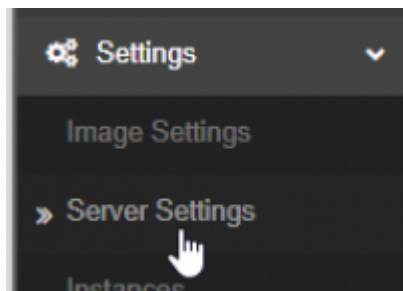

Email settings

Change settings about the email's sent by PMA.core

Setting	Description
Use Two Factor Authentication	Whether to enable two factor authentication for logging in
The name of the two factor authentication	The name of the keys generated for two factor authenticated
Email Address	The email present in all emails sent by PMA.core
Send emails using	Whether to sent emails using the SendGrid service or by a mail server
SendGrid Api Key	The api key to authenticate against SendGrid service

A screenshot of the 'Server Settings' web interface. It has a breadcrumb 'Home > Server settings' and a tabbed interface with 'Email Settings' selected. The form includes: a checkbox for 'Use Two Factor Authentication' (unchecked); a text field for 'The name of the two factor authentication' with value 'PMA Core'; a text field for 'Email Address' with value 'info@pathomation.com' and a note 'This address will show up as the sender of the emails'; a dropdown for 'Send emails using' with 'SendGrid' selected; and a text field for 'SendGrid Api Key' with a long alphanumeric string. At the bottom are buttons for 'Save', 'Clear settings', 'Test email', and 'Mail log'.

Mail log

A subsection of the [events page](#) showing only mail entries

Mail Log		
Time	Username	Description
5/25/2021, 11:59:09 AM	korinna	Mail sent with subject: PMA control - Registration Successful
6/7/2021, 2:01:03 PM	charlotte	Mail sent with subject: Password Reminder
7/12/2021, 12:36:02 PM	charlotte	Mail sent with subject: Password Reminder
8/9/2021, 9:58:10 AM	wordpress	Mail sent with subject: Password Reminder
2/16/2022, 3:14:50 PM	veerle	Mail sent with subject: Password Reminder
Showing 1 to 5 of 5 entries		
Previous 1 Next		

User settings

Settings regarding the user and password management

Setting	Description
Max invalid login attempts before lock out}	The max invalid login attempts allowed for a user. Use 0 to never lockout
Automatically lift lock out after (seconds)	Automatically lift lock out after (seconds)
Allow duplicate emails	Whether the system can allow users with the same email but different username
Force change password after days	Require users to change password after the specified amount of days
Password must be different from the last used	Require users to not use the same password twice

User Settings

Home > User settings

Email Settings

User settings

Password generation

LDAP Settings

OAuth Settings

Installation History

Max invalid login attempts before lock out

5

Automatically lift lock out after (seconds)

5

Allow duplicate emails

☒

Force change password after days

0

Use 0 to disable

Password must be different from the last used

0

Use 0 to disable

Save

Password generation

Settings about the automatic generator of passwords in the [Users -> Create](#) page

Setting	Description
Case	The type of the letters case
Length	The lenght of the generated password
Include letters	Whether to include letters
Include numbers	Whether to include numbers
Include special characters	Whether to include special characters like !@\$ etc

LDAP settings

Shows the list of configured LDAP servers used for authentication by PMA.core. You can also delete,update and create new LDAP server settings

Add new LDAP server

Clicking the *Add new LDAP server* button will take you to the page to fill the required info to add a new LDAP server that will be used for authentication by PMA.core.

Setting	Description
LDAP Server	The LDAP server url
LDAP distinguished name	The distinguish name of the LDAP directory to connect to
Username	The username for authentication against the LDAP server
Password	The password for authentication against the LDAP server

Add new LDAP Server

LDAP Server

LDAP distinguished name

Username

Password

Save

Back to List

OAuth settings

Shows the list of configured OAuth servers used for authentication by PMA.core. You can also delete,update and create new OAuth server settings

Setting	Description
OAuth callback url	The callback url required by OAuth authentication protocol by third party applications using PMA.core for authentication

OAuth Settings

OAuth Settings

Callback URL: https://demo.pathomation.com/pma-core/2.0.1/Layer/OAuthCallback

OAuthName	OAuth Client Id	OAuth Public Key	Auth url	Token url	Client	
Identity Server	sub		https://demo.pathomation.com/pma-core/2.0.1/Layer/OAuthCallback	https://demo.pathomation.com/pma-core/2.0.1/Layer/OAuthCallback	PMAcore	# Edit

Add new OAuth server

Add new OAuth server

Clicking the *Add new OAuth Server* button will take you to the page to fill the required info to add a new OAuth server that will be used for authentication by PMA.core. PMA.core uses the **Browser Flow** protocol for authentication as described the OAuth 2.0. For more information see [OAuth2.0](#)

Setting	Description
Name	The name of the OAuth server
Public Key	The public key for authenticating against OAuth server
User Id claim	The claim that contains the user id
Auth Url	The url used for Auth by the OAuth provider
Token Url	The url used for token access by the OAuth provider
Client Id	The client Id for PMA.core as configured in the OAuth provider

Setting	Description
Client Secret	The client secret for the corresponding client id

Add new OAuth Server

Name	<input type="text"/>	
Public Key	<input type="text"/>	The Public Key field is required.
User Id claim	<input type="text"/>	The User Id claim field is required.
Auth Url	<input type="text"/>	
Token Url	<input type="text"/>	
Client Id	<input type="text"/>	
Client Secret	<input type="text"/>	Leave empty if no client secret is required for code to token exchange
Callback Url	https://docs.pathomation.com/docs/pma.core/2.0.1/docs/OAuthCallback	
<input type="button" value="Save"/> <input type="button" value="Back to List"/>		

Installation history

Shows a brief history of version changes and updates to the PMA.core server

Installation History

Home > Installation History

Email Settings

User settings

Password generation

LDAP Settings

OAuth Settings

Installation History

Date	Version	API Version	PMA-UI Version	License	Download	Request code	License change date
3/26/2022 5:17:38 PM	2.0.2.2022	1.3.8	2.16.0	IS	IS	Download PMA-UI (2.16.0) (2022-03-26)	
2/13/2022 9:23:44 AM	2.0.1.2022	1.3.8	2.16.0	IS	IS	Download PMA-UI (2.16.0) (2022-02-13)	
2/13/2022 9:20:47 AM	2.0.1.1902	1.3.8	2.16.0	IS	IS	Download PMA-UI (2.16.0) (2022-02-13)	

From:

<https://docs.pathomation.com/pma.core/2.0.1/> - **PMA.core 2.x**

Permanent link:

https://docs.pathomation.com/pma.core/2.0.1/doku.php?id=settings_server&rev=1648824906

Last update: **2022/04/01 17:55**

