


# Context menu

The product indication label in the top-left corner of the PMA.studio and PMA.vue user interface serves as a drop-down menu. This is the starting tab of the ribbon and provides set of menu items used by the application.

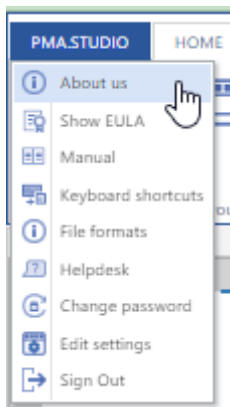

The items in the menu are the same for PMA.studio and PMA.vue.

## About us

This menu option brings up the about dialog with general information about the system.

This menu option is helpful in determining whether you use PMA.vue or PMA.studio, and what version.



If connected to a PMA.core instance and looking at a slide, the information is included in the dialog (as in the example above).

## Show EULA

Pulls up the End-User License Agreement.



An organization can customize the EULA based on specific needs or use cases (e.g. [RUO](#) related).

## Manual

Takes you to [this website](#).

## Keyboard shortcuts

Brings up available keyboard shortcuts for navigating between slides and within a slide.



Keyboard shortcuts that were defined in addition to make certain ribbon-options more directly accessible are not included in this list.

## File formats

This menu option takes you to [a live website where all currently supported file formats are listed](#).

### Supported WSI formats

- Number of brightfield formats: 39
- Number of fluorescent formats: 19
- Number of z-stacking formats: 19

| VENDOR / FILE FORMAT | EXTENSION   | COMPRESSION                            | BRIGHTFIELD | FLUORESCENCE | Z-STACKING |
|----------------------|-------------|--|-------------|--------------|------------|
| TIFF                 | .tif, .tiff | JPEG, JPEG2000, LZW, Deflate, Raw, RLE | ✓           |              |            |
| JPEG                 | .jpeg, .jpg | JPEG                                   | ✓           |              |            |
| JPEG 2000            | .jp2        | JPEG2000                               | ✓           |              |            |
| PNG                  | .png        | PNG                                    | ✓           |              |            |
| Olvminus VSI         | .vsi        | lossless JPEG, JPEG, ...               | ✓           | ✓            | ✓          |

The list may slightly differ from the actual file formats that *you* can read, in case you are not using the latest [PMA.core tile server](#) version.


## Helpdesk

Pathomation has a [helpdesk](#), where anybody with any question about us or our products, can file a ticket 24/7.

We then do our best to get back to you as soon as possible, with an adequate solution.

## Change password

If you're logged in to a particular PMA.core instance (PMA.vue users only see one anyway; PMA.studio may have multiple to choose from), you can use this dialog to change your password.

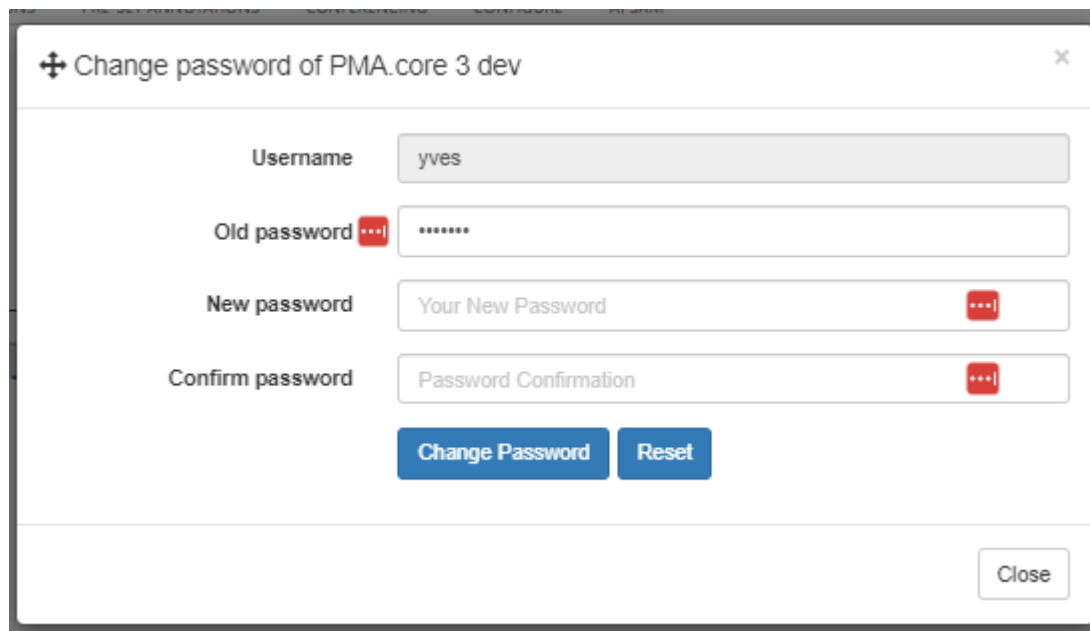


The image shows a dialog box titled "Change password of PMA.core 3 dev". It contains four input fields: "Username" with the value "yves", "Old password" (masked with dots), "New password" (with placeholder "Your New Password"), and "Confirm password" (with placeholder "Password Confirmation"). There are two buttons at the bottom: "Change Password" and "Reset". A "Close" button is located in the bottom right corner.

Passwords are case sensitive! If your password is Gr33nGr4ss, you may only enter it exactly so. Alternatives like gr33ngr4ss or GR33NGR4SS will not allow you access to the system.

## Edit settings

Edit settings specifically linked to your own profile.

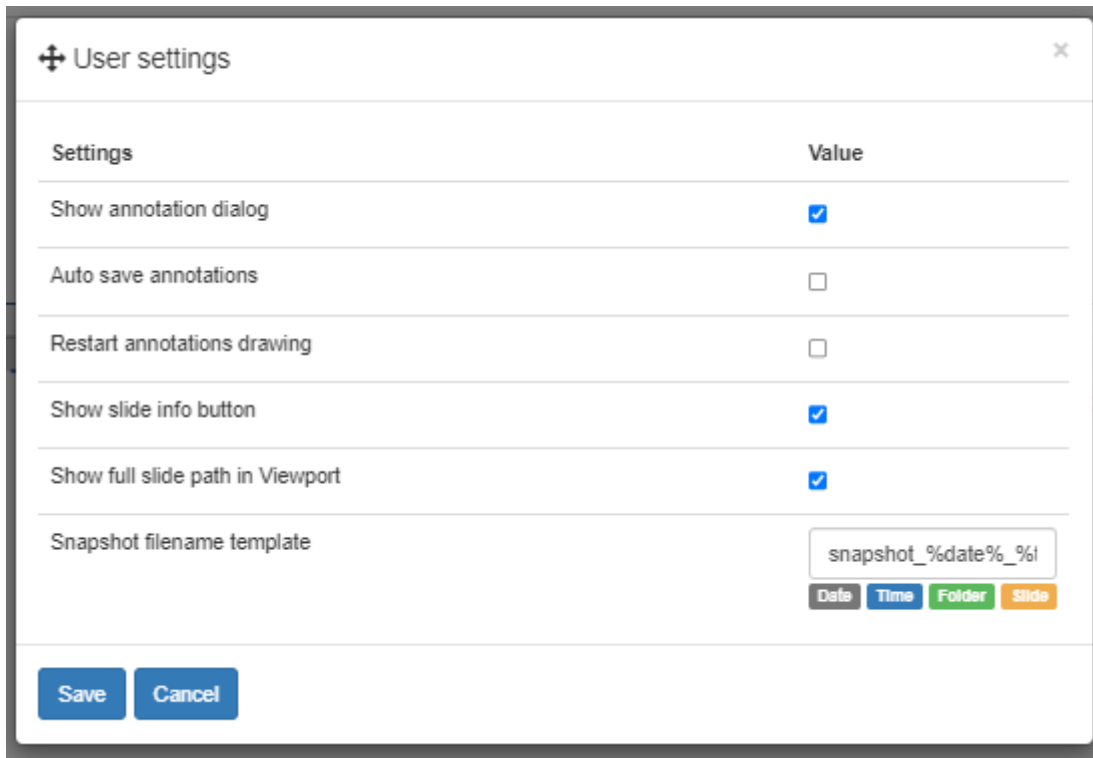

Note that your profile is linked to PMA.studio, and not to one particular PMA.core instance!

## Sign out

Sign out of all connected PMA.core instances and return to the login screen.

From:  
<https://docs.pathomation.com/pma.studio/2.1.0/> - **PMA.studio 2.1.0**

Permanent link:  
[https://docs.pathomation.com/pma.studio/2.1.0/doku.php?id=context\\_menu&rev=1676026222](https://docs.pathomation.com/pma.studio/2.1.0/doku.php?id=context_menu&rev=1676026222)

Last update: **2023/02/10 13:50**

