

# Working with grids

PMA.studio uses viewport as a grid system. UI provides a separate tab “GRID” to perform grid operations on the viewport.




“GRID” tab consists of 6 sections:

1. Size
2. Synchronize
3. R/C
4. Slide
5. Snapshot
6. Share


## Size

Size section allows you to break viewport as a grid system with different sizes. This includes grid sizes of 1x1,2x2,3x3 and even you can customized the grid size according to your need.

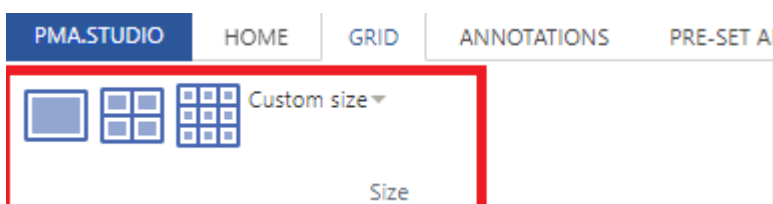


### (1x1 Grid size)

First grid option allows you to divide the viewport into 1x1 column grid.

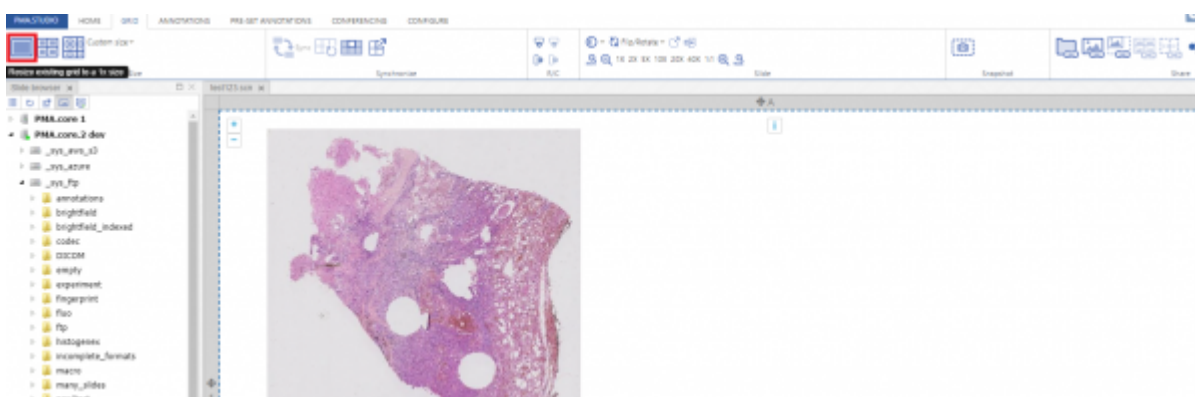

### (2x2 Grid size)

Second grid option divides the viewport into 2x2 column grid.

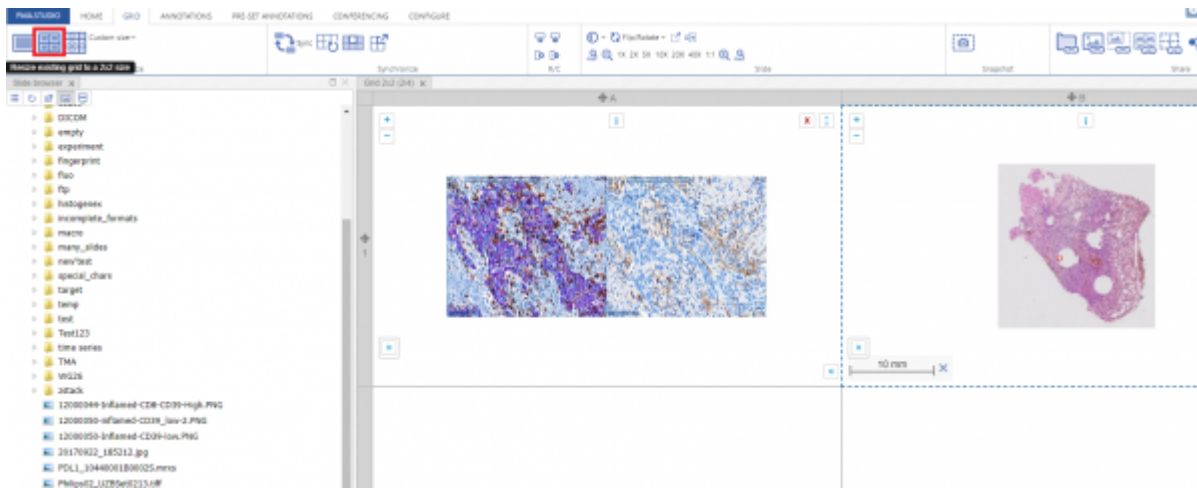

### (3x3 Grid size)

Third grid option divides the viewport into 3x3 column grid.



### (Custom Grid size)

You can even customize the grid according to your need by using the custom grid option.



### Synchronize

From:

<https://docs.pathomation.com/pma.studio/2.1.0/> - **PMA.studio 2.1.0**

Permanent link:

[https://docs.pathomation.com/pma.studio/2.1.0/doku.php?id=slides\\_grid&rev=1650995506](https://docs.pathomation.com/pma.studio/2.1.0/doku.php?id=slides_grid&rev=1650995506)

Last update: **2022/04/26 20:51**

