


Edit Training Session

In the Training Section page - accessible via the sidebar - you can view all existing training sessions [see setup a session for how to create one](#). Clicking on 'Edit' will allow you to edit one.

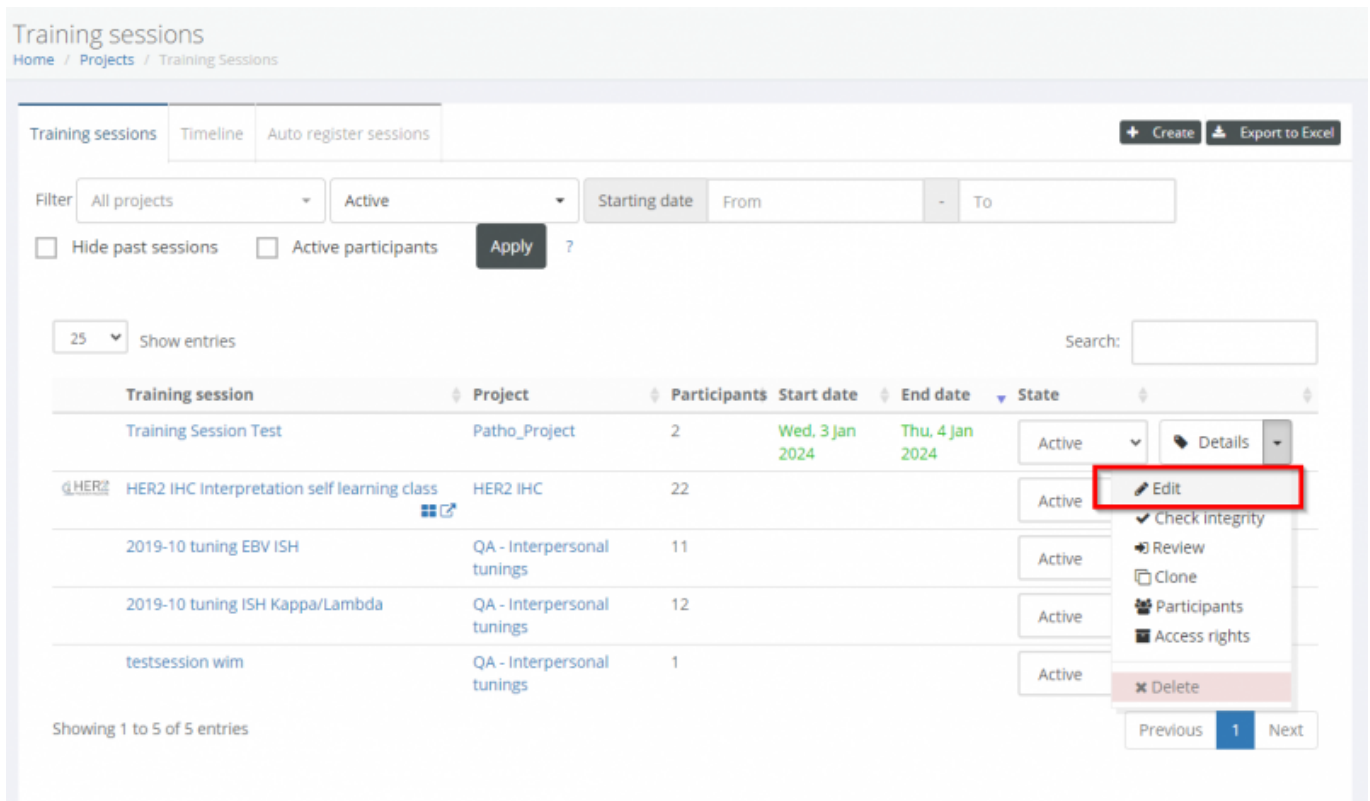

The screenshot displays the 'Training sessions' management page. At the top, there are tabs for 'Training sessions', 'Timeline', and 'Auto register sessions'. Below the tabs are filter options for 'All projects', 'Active', and 'Starting date' (From - To). There are also checkboxes for 'Hide past sessions' and 'Active participants', and an 'Apply' button. A search bar is located on the right. The main content is a table with columns: Training session, Project, Participants, Start date, End date, and State. The table lists five sessions, with the second one, 'HER2 IHC Interpretation self learning class', having an open context menu. The menu options are: Edit (highlighted with a red box), Check integrity, Review, Clone, Participants, Access rights, and Delete. The table footer shows 'Showing 1 to 5 of 5 entries' and pagination controls.

Training session	Project	Participants	Start date	End date	State
Training Session Test	Patho_Project	2	Wed, 3 Jan 2024	Thu, 4 Jan 2024	Active
HER2 HER2 IHC Interpretation self learning class	HER2 IHC	22			Active
2019-10 tuning EBV ISH	QA - Interpersonal tunings	11			Active
2019-10 tuning ISH Kappa/Lambda	QA - Interpersonal tunings	12			Active
testsession wim	QA - Interpersonal tunings	1			Active

Once in the edit menu you will be able to:

1. Insert an image
2. View the associated project
3. Change the name of the session
4. Check/uncheck if this is a course
5. Insert a short description of the training session
6. Define or edit the start and end date of the training session's availability
7. Write a welcome message for Users
8. Review / preview your changes

Editing training session Training Session Test

Home / Projects / Patho_Project / Training Sessions / Training Session Test

Details	Edit	Participants	Consensus data	Email	Integrity check	Slidebox	Certification	8 Review
---------	-------------	--------------	----------------	-------	-----------------	----------	---------------	---------------------------------

1	Logo	<input type="button" value="No file selected"/> Max file size: 10Mb		
2	Project	Patho_Project		
3	Title	<input type="text" value="Training Session Test"/>		
4	Course mode	<input type="checkbox"/>		
5	Short description	<div style="border: 1px solid #ccc; height: 100px;"></div>		
6	Starts on	<input type="text" value="2024-01-03"/>	Ends on	<input type="text" value="2024-01-04"/>
7	Welcome message	<div style="border: 1px solid #ccc; padding: 5px;"><p> 14 B <i>I</i> <u>U</u> S A </p><div style="border: 1px solid #ccc; height: 100px;"></div></div>		

Further options include:

9. Add case collection notes
10. Select a reference image collection
11. Select which case collections associated with the Project should be included/excluded from this training session with the checkboxes
12. Change what role types can share the training session via a dropdown menu
13. Select which (if any) questionnaire should be displayed to users prior to completing the case collection content
14. Select which server/geo-replication the session should be hosted from
15. Select whether users can see consensus data or not
16. Check/uncheck attributes that an admin has defined [see section on system attributes](#)

9 Case collection notes

10 Reference images collection

▼

11 Case collections ?

Test

Test 3

Test 2

12

Allow sharing

Managers & Supervisors ▼

13

Questionnaire form

Test form ▼

14

Region

PMA.core ▼

15

User can toggle consensus data

16

Attributes

Billable

Free text notes

From:

<https://docs.pathomation.com/pathotrainer/> - Pathotrainer

Permanent link:

https://docs.pathomation.com/pathotrainer/doku.php?id=edit_a_training_session&rev=1704294924

Last update: 2024/01/03 18:15

