


## Home: Folder Tools


The Folder Tools group in the Home tab contains controls for searching & navigating folders

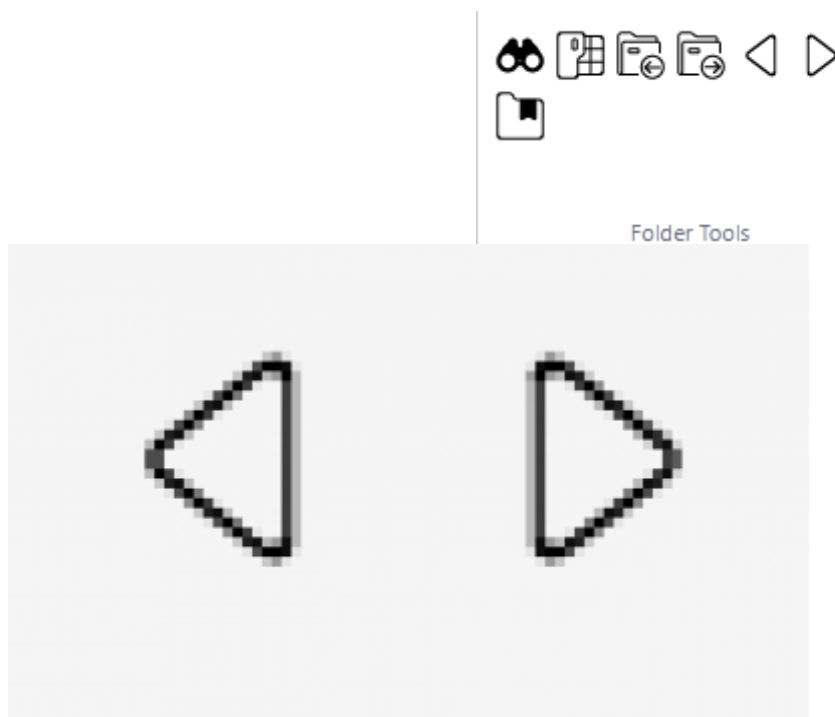

The icons, from left to right

### Top Row

1. - Search icon - [read more about searching here](#)
2. - Open current folder in grid - opens the menu below:

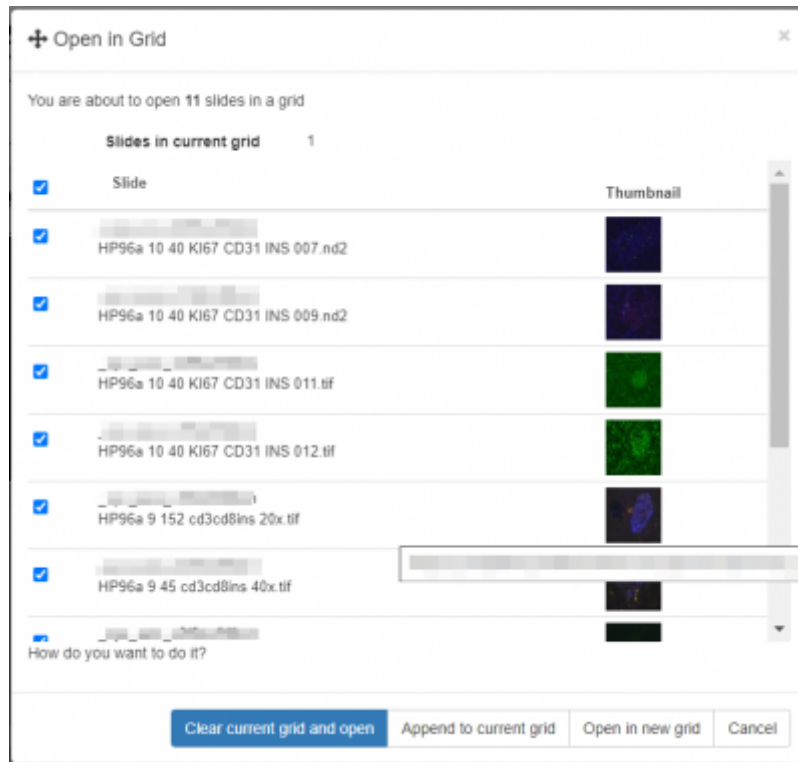

Here you can choose images from the currently highlighted folder in the slide browser to open into your current grid, or make a new grid.

- 3. - Previous folder - selects and opens the previous folder
- 4. - Next folder - selects and opens the next folder
- 5. - Previous file - opens previous image in folder 6. - Next file - opens next image in folder

From:  
<https://docs.pathomation.com/pimsls/3.0.1/> - **PIMS LS 3.0.1**

Permanent link:  
[https://docs.pathomation.com/pimsls/3.0.1/doku.php?id=folder\\_tools&rev=1746897758](https://docs.pathomation.com/pimsls/3.0.1/doku.php?id=folder_tools&rev=1746897758)

Last update: **2025/05/10 20:22**

