

Configure user accounts

Add new user accounts

To add a user click on the “Add” button on the top right of the table.



The following fields can be defined for each user:

A screenshot of the 'Add user' form. The form has a title 'Add user' and a breadcrumb 'Home > Users > Create'. It contains several input fields: 'First Name', 'Last Name', 'Login', 'Password' (with a 'Generate' button), 'Confirm password', 'Email', and 'Description'. At the bottom, there are three checkboxes: 'Administrator', 'Can Annotate', and 'Suspended'.

Field	Description
First Name	The user's first name
Last Name	The user's last name
Login	The username, it must be unique irrespective of its casing, if “UsEr” exists “uSeR” cannot exist
Password	The user's password
Confirm Password	Insert the password again for confirmation
Email	The user's email. Must be unique if the corresponding server setting is enabled
Administrator	Whether the user is an administrator, can login to the administrator UI and use the Admin API
Can Annotate	Whether the user can create new annotations on slides
Suspended	Whether the user is suspended, i.e. cannot authenticate and perform any action on the server

Modify user accounts

After creating a user you can click on the *Edit* button to modify an existing user. The page and fields are exactly the same as the *Create User* page (see previous paragraph). The only difference in this page is the ability **to leave the password/confirm password fields empty**, and the server will not change the user's password.

Other actions



Clicking on the down arrow in the user's list page you can see additional actions that can be performed on a user. Those actions are:


- Remove: This action will delete the user
- Remind password: Will send an email to the user with his password
- Reset password: This will set a new random generated password for this user and then sent an email containing the new password
- Unlock: Unlocks a locked out user, after too many unsuccessful login attempts ([see user states](#))

Audit trailing

Clicking on the Audit Trail tab of a user will display a chronological list of all changes performed to a user, as well as all changes performed to the access list of that user ([see access control lists](#))

This list contains the following useful info about a change in chronological order:

- ActionType: The type of a change i.e. whether the user was created, updated or deleted
- EventDate: The date the change occurred
- UserLogin: The user that performed the change
- List of user's field: The previous value for this user, and the changed value

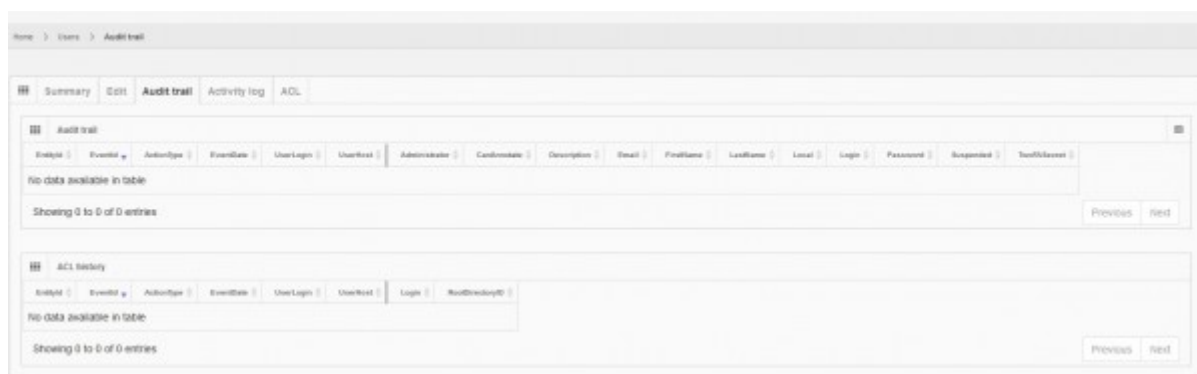

Activity Log

Clicking on the Activity Log tab will display a comprehensive list of all activities performed by the user over time

User log

Home > Users > Log

Summary

Edit

Audit trail

Activity log

ACL

System Access Log

Session ID	First Accessed	Last Accessed	Host Address	Host Name	Application	Success	Failure Reason	Source
h2VRoaAFVUKjgpyl6t04Zw2	04/04/2022 04:06:22	04/04/2022 16:06:16	85.73.193.76	85.73.193.76	PMA.core	true		Local

Showing 1 to 1 of 1 entries

Previous1Next

Slide Access Log

First Accessed	Last Accessed	Session ID	Slide Path	Slide Name
04/04/2022 15:06:26	04/04/2022 15:06:26	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refanncolations	CMJ-1 JP2K-03005 sys
04/04/2022 14:02:04	04/04/2022 14:02:04	h2VRoaAFVUKjgpyl6t04Zw2	_sys_sys_x3/test_js	appendicitis sys
04/04/2022 12:05:56	04/04/2022 12:05:56	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refbnglthed/Hanamabuu	2020-11-18 15:17:32 edpi
04/04/2022 12:03:36	04/04/2022 12:03:36	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refbnglthed/Hanamabuu	CMJ-1 ndpi
04/04/2022 12:00:19	04/04/2022 12:00:19	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refbnglthed/Hanamabuu	uric_uricfr_istatno2_115 Bone Marrow Gamma ndpi
04/04/2022 11:59:44	04/04/2022 11:59:44	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refbnglthed/Hanamabuu	H B2 Breast ndpi
04/04/2022 11:59:27	04/04/2022 11:59:27	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refbnglthed/Apeno (Lutka)	CMJ-1 sys
04/04/2022 11:57:11	04/04/2022 11:57:11	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refbnglthed/CHstachn2 2	B-1793193 HE annotations
04/04/2022 11:36:30	04/04/2022 11:36:30	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refbnglthed/CHstachn2 2	B-1793193 HE rnsa
04/04/2022 09:47:36	04/04/2022 13:50:26	h2VRoaAFVUKjgpyl6t04Zw2	_sys_refanncolations	CMJ-1 ndpi

Showing 1 to 10 of 13 entries

Previous12Next

From:
<https://docs.pathomation.com/pma.core/2.0.0/> - **PMA.core 2.x**

Permanent link:
https://docs.pathomation.com/pma.core/2.0.0/doku.php?id=user_config&rev=1649086259

Last update: **2022/04/04 18:31**

