

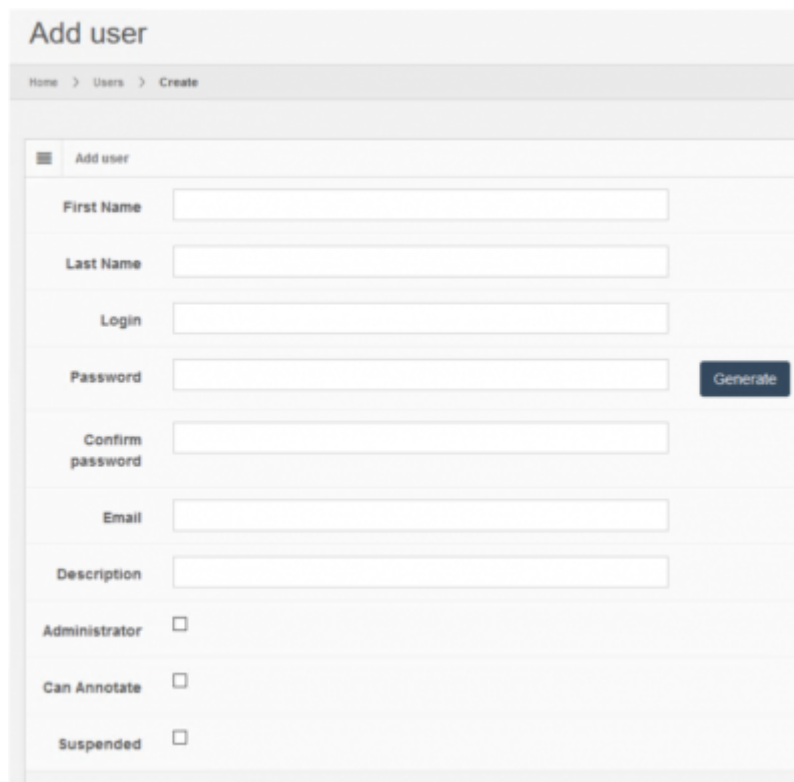
# Configure user accounts

## Add new user accounts

To add a user click on the “Add” button on the top right of the table.



The following fields can be defined for each user:

A screenshot of the 'Add user' form in a web application. The form is titled 'Add user' and has a breadcrumb trail 'Home > Users > Create'. It contains several input fields: 'First Name', 'Last Name', 'Login', 'Password' (with a 'Generate' button next to it), 'Confirm password', 'Email', and 'Description'. At the bottom, there are three checkboxes: 'Administrator', 'Can Annotate', and 'Suspended'.

Field	Description
First Name	The user's first name
Last Name	The user's last name
Login	The user's username. The username must be unique, irrespective of its casing. If user "UsEr" is already defined, you cannot also add user "uSeR".
Password	The user's password
Email	The user's email address
Administrator	Boolean value that indicates whether or not the user can login to the PMA.core administration interface
Can Annotate	Boolean value that indicates whether or not the user can create annotations for virtual slides
Suspended	Boolean value that indicates whether or not the user can login to the PMA.core and PMA.view systems

## Modify user accounts

## Delete user accounts

## Audit trailing

From:

<https://docs.pathomation.com/pma.core/2.0.1/> - **PMA.core 2.x**

Permanent link:

[https://docs.pathomation.com/pma.core/2.0.1/doku.php?id=user\\_config&rev=1646673092](https://docs.pathomation.com/pma.core/2.0.1/doku.php?id=user_config&rev=1646673092)

Last update: **2022/03/07 20:11**

