


Active sessions

In the session management section the currently active users are displayed. The system gives the ability to terminate a user's session.

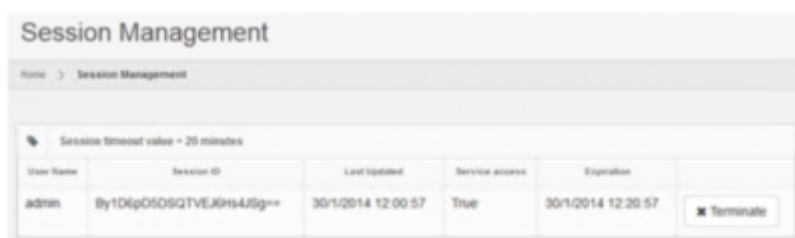

The screenshot shows a web interface titled "Session Management". Below the title is a breadcrumb "Home > Session Management". A message states "Session timeout value = 20 minutes". Below this is a table with the following columns: "User Name", "Session ID", "Last updated", "Service access", and "Expiration". There is a "Terminate" button next to the table. The table contains one row for a user named "admin" with session ID "8y1D6pO5OGQTVEJ6Hs4Z8g==", last updated "30/1/2014 12:00:57", service access "True", and expiration "30/1/2014 12:20:57".

User Name	Session ID	Last updated	Service access	Expiration	
admin	8y1D6pO5OGQTVEJ6Hs4Z8g==	30/1/2014 12:00:57	True	30/1/2014 12:20:57	Terminate

The lifespan of a session is governed by the ASP.Net session time out value. This value is displayed in the installation check page. When a session is idle for more time than what is defined as the session time out value, it becomes automatically invalid and the corresponding user will need to log in again.

You may delete an active session (e.g. when it's been idle for a large period of time; perhaps the end-user forgot to log out properly), but you cannot delete your own session.

From:

<https://docs.pathomation.com/pma.core/2.0.2/> - **PMA.core 2.x**

Permanent link:

https://docs.pathomation.com/pma.core/2.0.2/doku.php?id=active_sessions&rev=1591181307

Last update: **2020/06/03 13:48**

