


## Active sessions

In the session management section the currently active users are displayed. The system gives the ability to terminate a user's session.

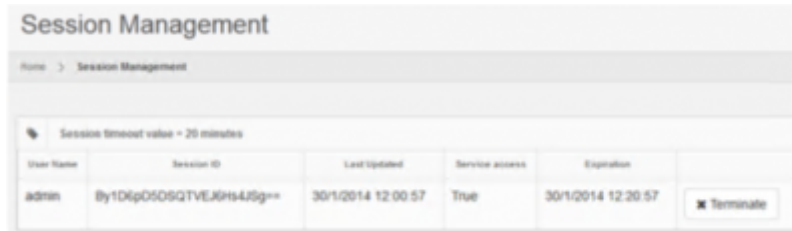

The screenshot shows a web interface titled "Session Management". Below the title, there is a breadcrumb "Home > Session Management". A message indicates "Session timeout value = 20 minutes". A table lists active sessions with columns for User Name, Session ID, Last updated, Service access, and Expiration. A "Terminate" button is visible next to the session for user "admin".

User Name	Session ID	Last updated	Service access	Expiration	
admin	By1D6p0S09QTVEJHs4J9g==	30/1/2014 12:00:57	True	30/1/2014 12:20:57	<input type="button" value="Terminate"/>

The lifespan of a session is governed by the ASP.Net session time out value. This value is displayed in the installation check page. When a session is idle for more time than what is defined as the session time out value, it becomes automatically invalid and the corresponding user will need to log in again.

You may delete an active session (e.g. when it's been idle for a large period of time; perhaps the end-user forgot to log out properly), but you cannot delete your own session.

From:  
<https://docs.pathomation.com/pma.core/2.0.2/> - **PMA.core 2.x**

Permanent link:  
[https://docs.pathomation.com/pma.core/2.0.2/doku.php?id=active\\_sessions&rev=1646664871](https://docs.pathomation.com/pma.core/2.0.2/doku.php?id=active_sessions&rev=1646664871)

Last update: **2022/03/07 17:54**

