


User monitoring

Overview table

User details view

By clicking on the user's login in the overview table, you can see a summary of the user's settings and activity. This page provides an overview of the user's information and the root directories they can access.

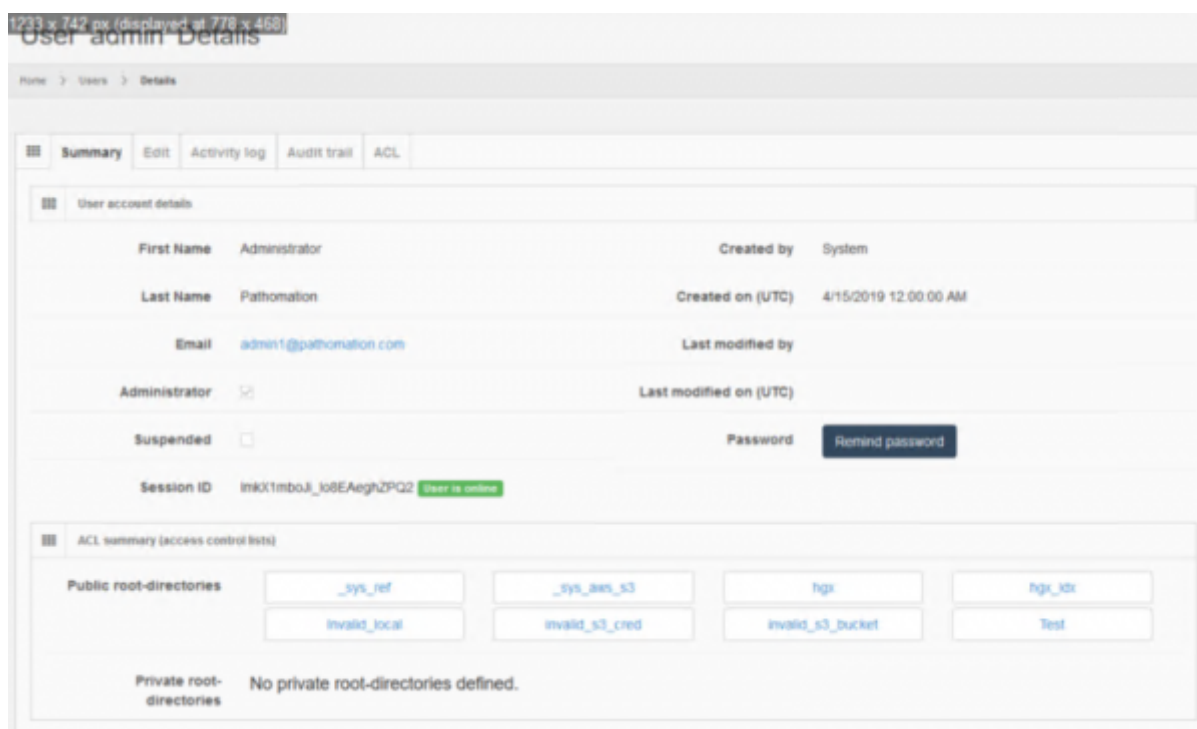

The screenshot displays the 'User admin Details' page. At the top, there's a breadcrumb trail: 'Home > Users > Details'. Below this is a tabbed interface with 'Summary', 'Edit', 'Activity log', 'Audit trail', and 'ACL'. The 'Summary' tab is active, showing 'User account details'. The details include: First Name (Administrator), Last Name (Pathomation), Email (admin1@pathomation.com), Administrator status (checked), Suspended status (unchecked), Session ID (lmkX1mboJl_lo8EAeghZPQ2), Created by (System), Created on (UTC) (4/15/2019 12:00:00 AM), Last modified by, and Last modified on (UTC). A 'Remind password' button is visible. Below the account details is the 'ACL summary (access control lists)' section. It shows 'Public root-directories' with a list of directories: _sys_ref, _sys_aes_s3, hgi, hgi_xdr, invalid_local, invalid_s3_cred, invalid_s3_bucket, and Test. The 'Private root-directories' section states 'No private root-directories defined.'

Access rights

Activity log

The activity log section provides a detailed list of events that involve the displayed user.

1227 x 817 px (displayed at 772 x 514)

User admin Log

Home > Users > Log

Summary Edit Activity log Audit trail ACL

System Access Log

Session ID	First Accessed	Last Accessed	Host Address	Host Name	Application	Success	Failure Reason
ImkX1mboJl_lo8EAeghZPQ2	4/15/2019, 9:50:49 AM	4/15/2019, 2:15:31 PM	85.73.248.94	85.73.248.94	PMA core	true	
	4/15/2019, 4:34:35 AM	4/15/2019, 4:34:35 AM	127.0.0.1	127.0.0.1	test	false	Invalid username or password
upjGsSqyVKQpikpbK268IA2	4/15/2019, 4:32:52 AM	4/15/2019, 9:45:22 AM	85.73.248.94	85.73.248.94	PMA core	true	

Showing 1 to 3 of 3 entries

Previous1Next

Slide Access Log

First Accessed	Last Accessed	Session ID	Slide Path	Slide Name
4/15/2019, 12:49:17 PM	4/15/2019, 2:15:32 PM	ImkX1mboJl_lo8EAeghZPQ2	._sys_ref/brightfield/Smartzoom	CMU-1.sz
4/15/2019, 12:48:25 PM	4/15/2019, 2:15:32 PM	ImkX1mboJl_lo8EAeghZPQ2	._sys_aws_s3/brightfield/Hanamatsu/CMU-1-VMS	CMU-1-40x - 2010-01-12 13:24:05.vms
4/15/2019, 12:47:59 PM	4/15/2019, 2:15:32 PM	ImkX1mboJl_lo8EAeghZPQ2	._sys_ref/brightfield/Hanamatsu/CMU-1-VMS	CMU-1-40x - 2010-01-12 13:24:05.vms
4/15/2019, 12:47:57 PM	4/15/2019, 2:15:32 PM	ImkX1mboJl_lo8EAeghZPQ2	._sys_aws_s3/brightfield/Generic_Png	png01_julia_4K.png
4/15/2019, 12:47:54 PM	4/15/2019, 2:15:32 PM	ImkX1mboJl_lo8EAeghZPQ2	._sys_ref/brightfield/Generic_Png	png01_julia_4K.png

Audit trail

From:
<https://docs.pathomation.com/pma.core/2.0.2/> - **PMA.core 2.x**

Permanent link:
https://docs.pathomation.com/pma.core/2.0.2/doku.php?id=user_monitoring&rev=1646672261

Last update: **2022/03/07 19:57**

