

User monitoring

Overview table


- The “Last login” column displays each user's last successful login date.
- The “Local” column indicates whether or not a user is managed by PMA.core or LDAP.
- The “Administrator” column indicates whether or not a user can access PMA.core. Non administrator users cannot perform system management tasks within PMA.core.
- The “Can Annotate” column indicates whether or not a user is allowed to create and edit annotations.
- The “Suspended” column indicates whether or not a user is currently allowed to login to either PMA.core or downstream applications. A suspended login is treated as invalid by the system.
- The “Locked out” column indicates whether or not a user is currently locked out due to consecutive invalid log in attempts.
- Clicking on a user's email address launches your default email client.
- Clicking on a user's login name displays a summary of the user's activity.

Next to each entry the following buttons appear:

Button	Description
Edit	Edits the selected user. This option allows the modification of the user's password, email and name. Furthermore the access level of the user can be changed (administrator or not) as well as whether or not the user is allowed to create and manage annotations. Finally a user can be (un)suspended within this option.
Delete	Deletes the selected user
Remind password	Sends a notification email to the selected user that contains their password in plain text
Reset and remind password	Sets a random password for the selected user and sends them a notification email that contains the new password in plain text
Unlock	Lifts a lock out.


User details view

By clicking on the user's login in the overview table, you can see a summary of the user's settings and activity. This page provides an overview of the user's information and the root directories they can access.

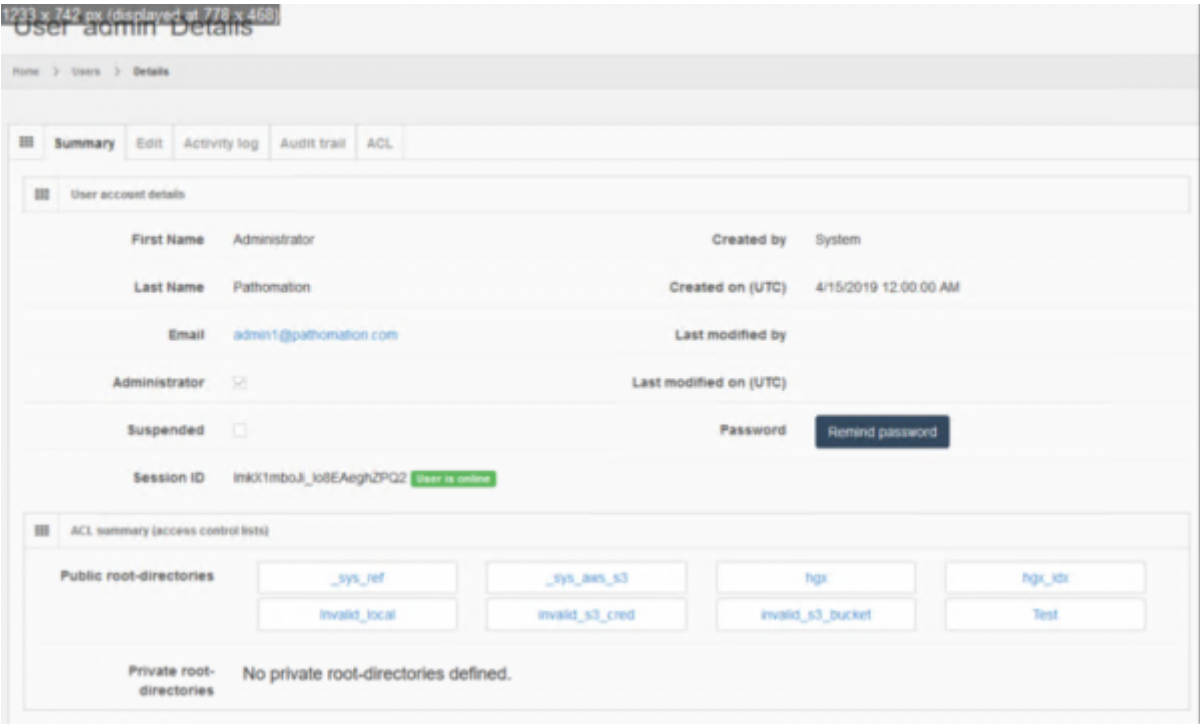

Access rights

The Access rights (ACL) section allows control over the [private root directories](#) the particular user may access. If no private root directories are defined, this list will be empty.

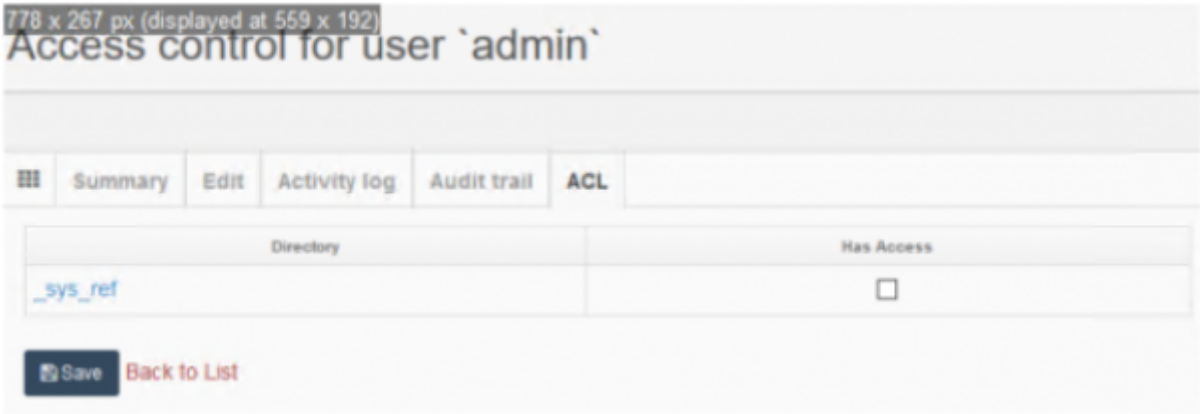

Activity log

The activity log section provides a detailed list of events that involve the displayed user.

1227 x 817 px (displayed at 772 x 514)

User admin Log

Home > Users > Log

SummaryEditActivity logAudit trailACL

System Access Log

Session ID	First Accessed	Last Accessed	Host Address	Host Name	Application	Success	Failure Reason
lnkX1mboJl_lo8EAeghZPQ2	4/15/2019, 9:50:49 AM	4/15/2019, 2:15:31 PM	85.73.248.94	85.73.248.94	PMA.core	true	
	4/15/2019, 4:34:35 AM	4/15/2019, 4:34:35 AM	127.0.0.1	127.0.0.1	test	false	Invalid username or password
upjGsSqVkiQpkpH268IA2	4/15/2019, 4:32:52 AM	4/15/2019, 9:45:22 AM	85.73.248.94	85.73.248.94	PMA.core	true	

Showing 1 to 3 of 3 entries

Previous1Next

Slide Access Log

First Accessed	Last Accessed	Session ID	Slide Path	Slide Name
4/15/2019, 12:49:17 PM	4/15/2019, 2:15:32 PM	lnkX1mboJl_lo8EAeghZPQ2	._sys_ref/brightfield/Smartzoom	cmu-1.sz
4/15/2019, 12:48:25 PM	4/15/2019, 2:15:32 PM	lnkX1mboJl_lo8EAeghZPQ2	._sys_aws_s3/brightfield/Hamamatsu/CMU-1-VMS	CMU-1-40x - 2010-01-12 13:24:05.vms
4/15/2019, 12:47:59 PM	4/15/2019, 2:15:32 PM	lnkX1mboJl_lo8EAeghZPQ2	._sys_ref/brightfield/Hamamatsu/CMU-1-VMS	CMU-1-40x - 2010-01-12 13:24:05.vms
4/15/2019, 12:47:57 PM	4/15/2019, 2:15:32 PM	lnkX1mboJl_lo8EAeghZPQ2	._sys_aws_s3/brightfield/Generic_Png	png01_julia_4K.png
4/15/2019, 12:47:54 PM	4/15/2019, 2:15:32 PM	lnkX1mboJl_lo8EAeghZPQ2	._sys_ref/brightfield/Generic_Png	png01_julia_4K.png

Audit trail

The audit trail section displays a log of modification events that involve the current user.

1230 x 471 px (displayed at 854 x 327)

Audit trail for user: test

Home > Users > Audit trail

SummaryEditActivity logAudit trailACL

Audit trail

EntryId	EventId	ActionType	EventDate	UserLogin	UserHost	Administrator	CardAnnotate	Description	Email	Firstname	Lastname
0	174	Update	4/15/2019, 4:34:15 PM	admin	85.73.248.94	False	True	Test user	test123@pathomation.com	Test	test
0	173	Create	4/15/2019, 4:34:00 PM	admin	85.73.248.94	False	False	Test user	test123@pathomation.com	Test	test

Showing 1 to 2 of 2 entries

Previous1Next

From:
<https://docs.pathomation.com/pma.core/2.0.2/> - PMA.core 2.x

Permanent link:
https://docs.pathomation.com/pma.core/2.0.2/doku.php?id=user_monitoring&rev=1646673115

Last update: 2022/03/07 20:11

