

Annotations

Basic Controls

To start with we have basic set of controls that allows us perform different operations on an annotation. It includes saving, toggling and styling of the annotation.



Toggling Annotation Panel

The very first button in the annotation panel is “Toggle Annotation Panel” which adds the annotation panel to the viewport. The panel shows details about all the annotations you added to a particular slide with different stats as well.



Saving and auto-saving annotations

Next is “Save Annotations” button that allows you to save annotation to the PMA.core server. Once you click on it you will get a notification if annotation was successfully save or not.




Controlling the visibility of annotations elements

(Toggle Annotations)


To control the visibility of annotation elements we can click on the “Toggle Annotations” button. The

settings on the ribbon are applied global to all annotations. You can control the visibility of individual annotations via the annotations panel.

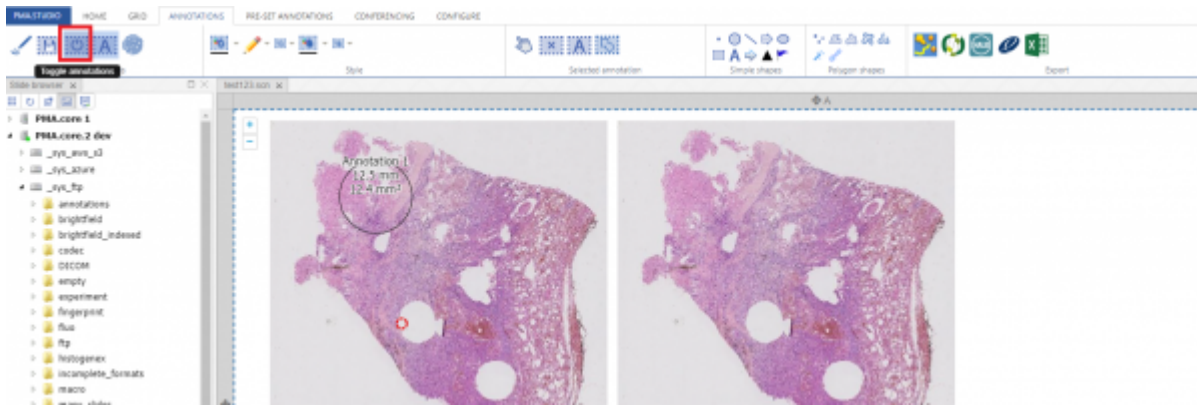


(Toggle Annotation Text)

“Toggle Annotation Text” button allows us to control visibility of annotation text only.

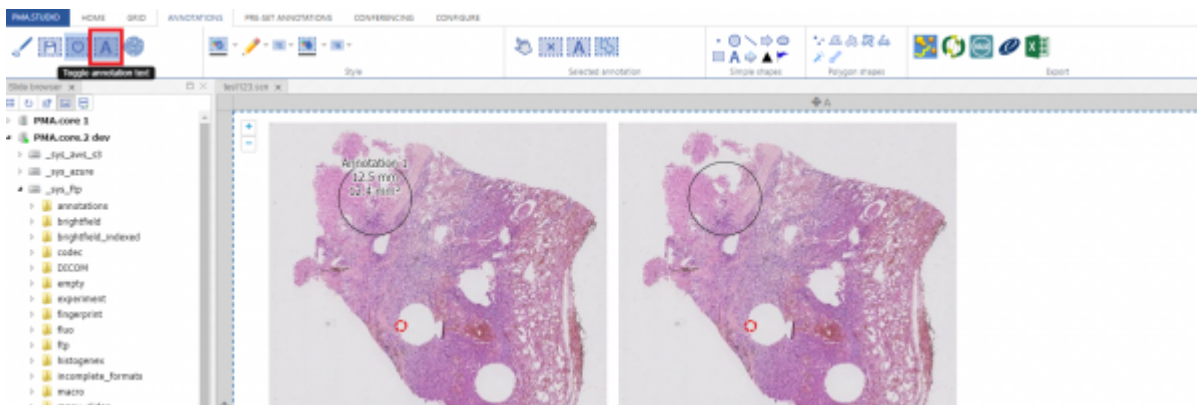

Styling your annotations

PMA.studio allows adding different styling to the annotations. You can change annotation color, pen width and opacity using these controls:

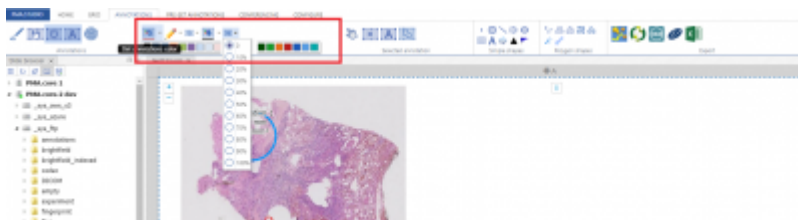

From:
<https://docs.pathomation.com/pma.studio/2.1.0/> - **PMA.studio 2.1.0**

Permanent link:
https://docs.pathomation.com/pma.studio/2.1.0/doku.php?id=annotations_configure&rev=1650553028

Last update: **2022/04/21 17:57**

