

# Edit annotations



Selected annotation

PMA.studio allows you to edit, delete individual, delete all, modify and merge annotations. Once you create annotations they are added to annotation tray panel which allows you to select an annotation and perform certain actions on it.

## Cursor




This button allows you to cancel any ongoing annotation operation. If you're in continuous drawing mode e.g. and want to stop drawing, you can click on this cursor icon and your mouse will go back to its default operational "pan and zoom" mode.


## Delete




This button delete the currently selected annotation. Once you click the button it shows you a confirmation popup and once you press "Delete" button, annotation gets deleted.

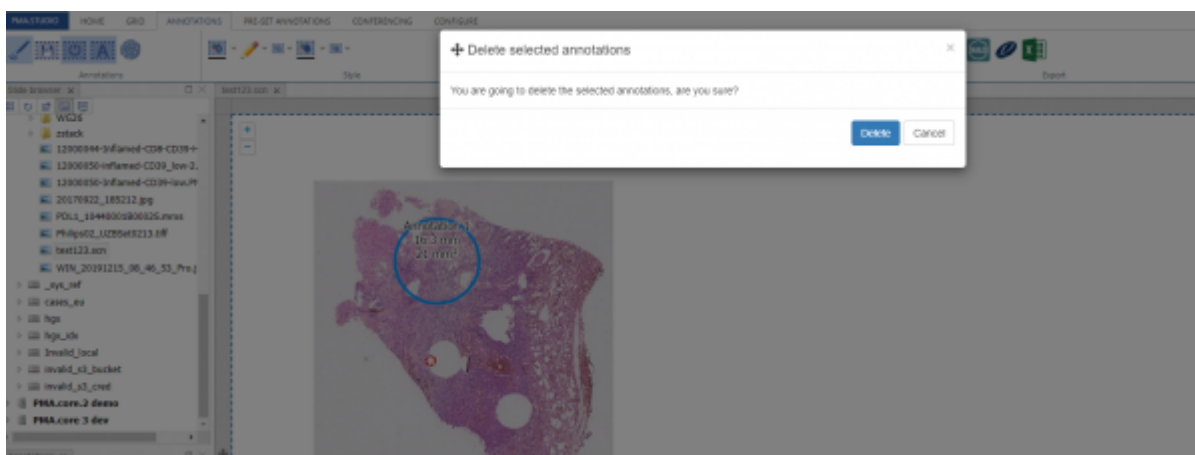

## Delete All Annotations

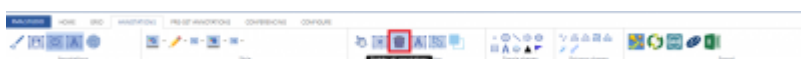


Delete all annotations button allows user to delete all annotations at once. Once you click the button it shows you a confirmation popup and once you press “Delete” button, all annotations get deleted.

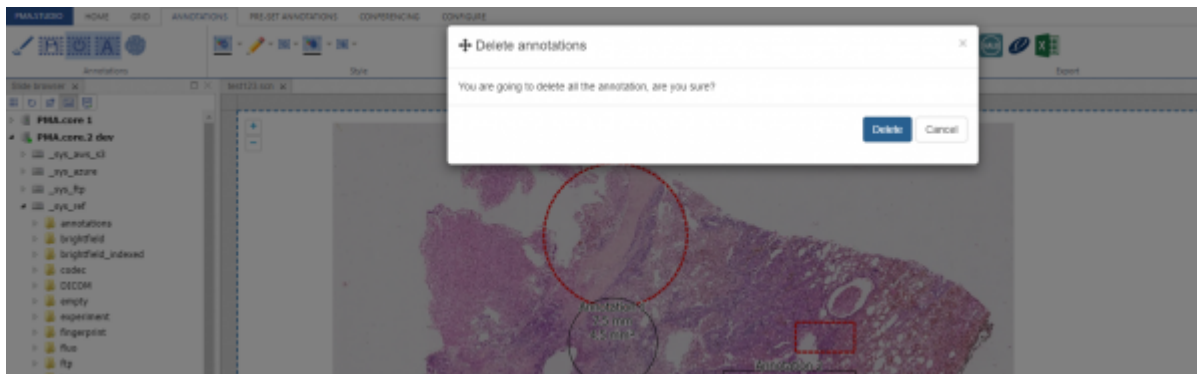



## Edit

“Edit” button allows you to edit an existing annotation. Once you select an annotation from annotation tray panel you see a popup window that allows you to change caption and classification of annotation. Once you modify the content and press “OK” annotation gets updated.

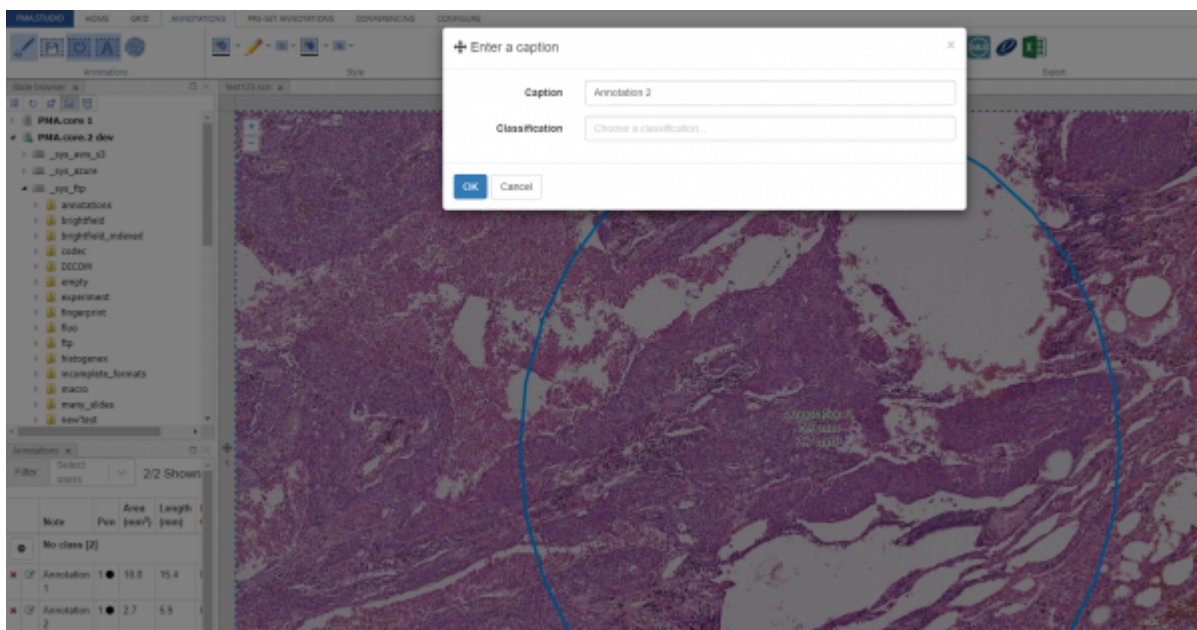


## Merge

Merge selected annotation button allows you to merge two annotations together. Right now PMA.Studio allows merging of “Line” and “Freehand” annotation types. To merge annotations you first select the two annotations using shift key and the press “merge selected” button and both annotations will be merged.

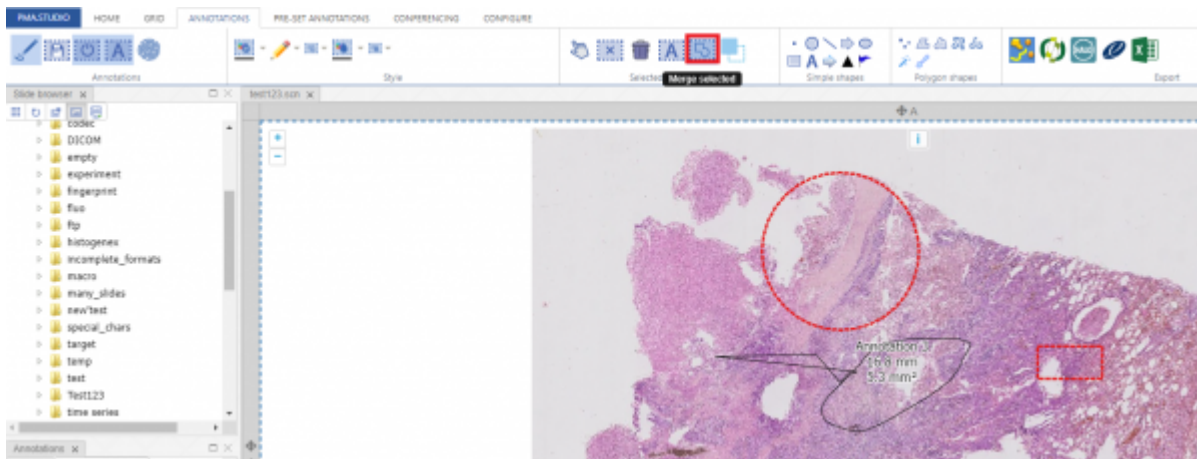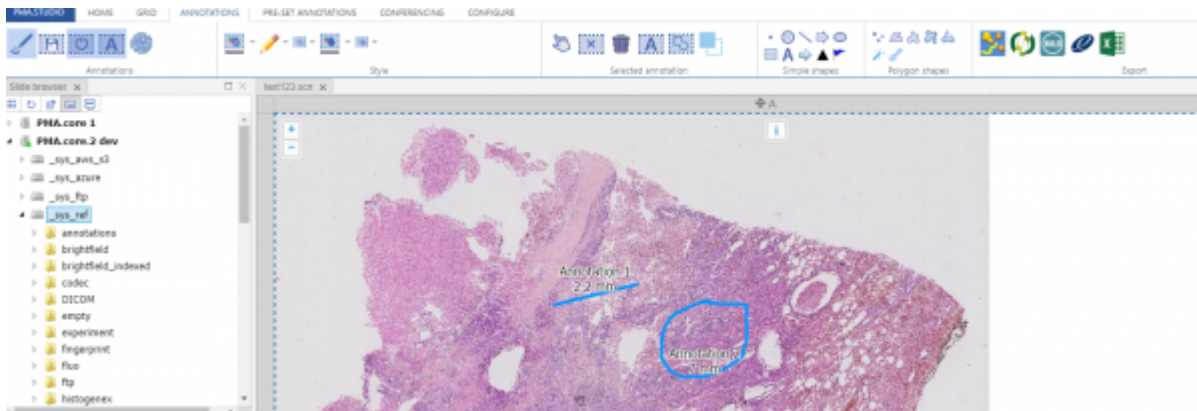

## Modify

Modify button allows user to modify an existing annotation. To modify an annotation you need to select the annotation and then press the “modify” button. Once you are done with the modification, you can save the annotation.

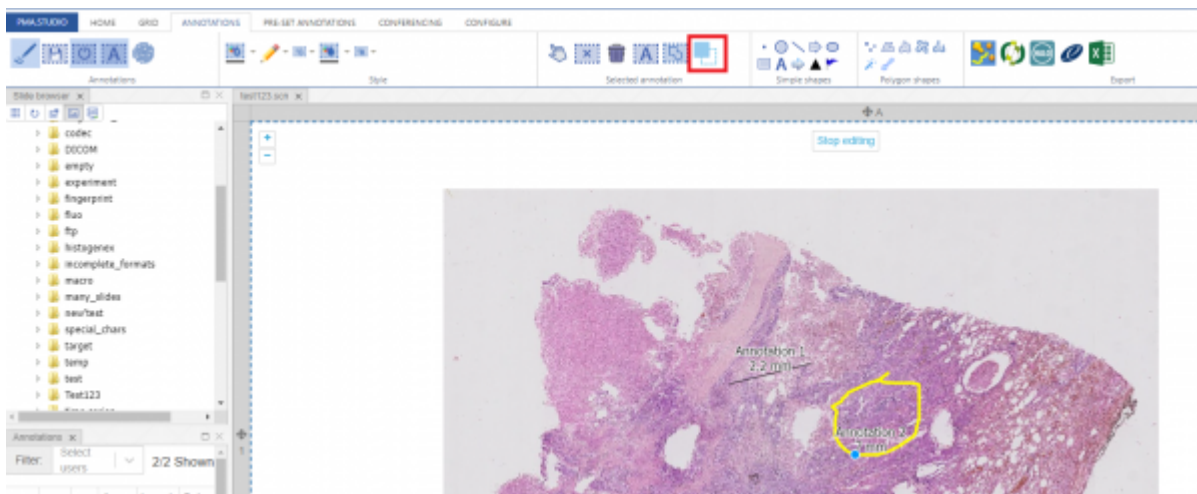

From:

<https://docs.pathomation.com/pma.studio/2.1.0/> - **PMA.studio 2.1.0**

Permanent link:

[https://docs.pathomation.com/pma.studio/2.1.0/doku.php?id=annotations\\_edit&rev=1672135008](https://docs.pathomation.com/pma.studio/2.1.0/doku.php?id=annotations_edit&rev=1672135008)

Last update: **2022/12/27 12:56**

